



**About Us** 



#### **Surgsci Medical Ltd.**

Helping medical professionals do their work is our work. Explore our range of medical solutions, design to improve performance at every point of care. Surgsci Medical Ltd. was established in Shenzhen. We have built many of our labs together with many institutes in China including medical polymer material laboratory, metal material laboratory, microbiology laboratory and surgical energy laboratory. We have an 80 people team with plentiful technical clinical experience. Since establishment, we have obtained more than 50 patents both domestic and overseas. In 2019, Surgsci Medical Ltd. was recognized as a national high-tech enterprise.





Click for more details

Introduction

 $\mathsf{Endopass}^{\scriptscriptstyle\mathsf{TM}}$ 

Quickpass™

Easypass™&Bluntpass™

Hassonpass™

Without Valve™

Trocar Model-Single set

**Trocar Model-Kit** 

More






#### Olive-Shaped Design

- Olive-shaped top better fits to the palm.
- Round top of the obturator makes the contact with palm more comfortable.

#### Detachable Switch Cap

- Multi-disc valve in detachable cap flexibly enfold the inserted instrument and achieve better sealing.
- Sealing valve are multiple folded with rich malleability. Keep good sealing while swinging freely with the inserted instrument.



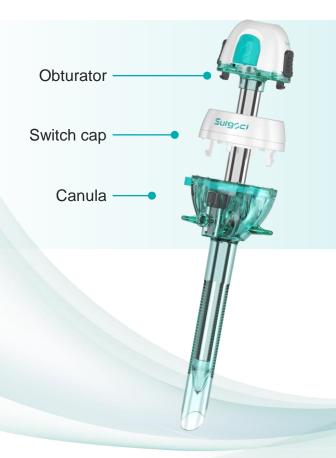

#### Fully Detachable structure

- Adopts fully detachable structure and composed of 3 parts with clipon link for assembly.
- Easy detached by pushing the button/handle and convenient for adjustment in Surgery.

#### Demo for fully detach

- Press button at both sides of obturator to release it.
- Stir the handle around canula body to separate the switch cap and cannula.
  Wider portal appears for specimen taken out during minimally surgery.







- Endopass visual laparoscopic trocar offers 7 sizes that can be used in various types of laparoscopic surgeries.
- Under the complexity of abdominal conditions, Endopass visual laparoscopic trocar gives visual feedbacks when puncturing with the help of endoscope to avoid the potential damage to the organs by mis-operation.



#### Funnel Lead-In

 Funnel shape Lead-In enabling smooth insertion, one-hand surgical instrument exchanging



#### Stable Thread Design

 Unique thread-shaped sleeve surface with inverted tooth design to increase sleeve fixation on abdominal wall and avoid sliding when the instruments entering or withdrawing from the sleeve.



#### None-bladed tip

Non-blade tip separates tissues instead of cutting.
Reduces the damage to abdominal wall and quickens the recovery after surgery



#### Endoscope secure lock


 Easy operating tunnel lock to secure the endoscope from sliding in order to provide stable visualization.





#### Available in 7 sizes

Code	Diameter	Length
T05075V	Ф5mm	75mm
T05100V	Ф5mm	100mm
T05150V	Ф5mm	150mm
T10100V	Ф10mm	100mm
T10150V	Ф10тт	150mm
T12100V	Ф12mm	100mm
T12150V	Ф12mm	150mm

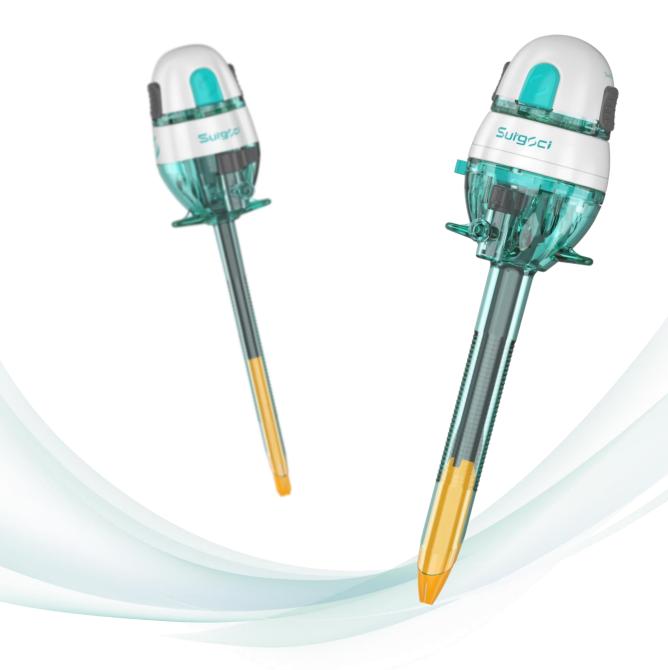

- Quickpass blade trocar provides 4 sizes that can be used in various types of laparoscopic surgeries.
- Automatically blade retracting system to ensure the retraction of the blade before and after insertion.
  Protecting both patient and surgeon from injuries.



#### Olive-shaped ergonomic design

 Olive shaped outline and Ergonomic design fits user' palm better. Make it more comfortable to hold and puncture.



#### Secure button

- Secure button released=Safe mode. Blade detained in sheath.
- Secure button pressed=Puncture mode. Blade pushed out when puncturing and retracted into sheath afterward to prevent injuries to organs.



#### Blade protective sheath design

 Unique blade design with sheath on tip. Preventing blade from abrasion and accidentally injury before insertion. Blade retracts back into sheath once entering peritoneum for safety operation.





#### Available in 4 sizes

Code	Diameter	Length
T05075P	Ф5mm	75mm
T05100P	Ф5mm	100mm
T10100P	Ф10mm	100mm
T12100P	Ф12mm	100mm



# Bluntpass™



- Easypass trocar provides 3 sizes and Bluntpass trocar provides 2.
- Non-blade puncture cone head separates tissue instead of cutting, reduces penetration damage and quickens recovery after surgery.



#### Finger grip

 Ergonomically designed finger grip fits user better. Allowing for easy rotating during insertion and offering a firm holding position.



#### Double seal structure

 Double sealing units design. Detachable multi-disc sealing structure and one-way duckbill sealing valve allows the entry of different sizes of instruments and maintain the pneumoperitoneum while operating



#### Easypass puncturing tip

 With plastic sharp edge on both sides of the tip, the unique designed obturator divides tissue easily and precisely, reduces the damage to abdominal wall during puncturing.



#### Bluntpass puncturing tip

 With blunt and mellow puncturing tip, it avoids damage to abdominal wall and suitable for second puncturing.



# Surgoci

#### Available in 3 sizes-Easypass trocar

Code	Diameter	Length
T05100S	Ф5mm	100mm
T10100S	Ф10тт	100mm
T12100S	Ф12mm	100mm

#### Available in 2 sizes-Bluntpass trocar

Code	Diameter	Length
T10100D	Ф10mm	100mm
T12100D	Ф12mm	100mm



- Hassonpass trocar provides 2 sizes for puncturing surgery
- The blunt and mellow obturator tip pushes the tissue softly and avoids damage to abdominal wall when puncturing.
- Adjustable suture pile on canula tube. Suture line twined and fastened to suture piles making the canula settled on abdominal wall.





#### Suture line pile

Adjustable suture pile on canula tube.
Suture line twined and fastened to suture piles making the canula settled on abdominal wall.



#### Blunt puncturing tip

 With blunt and mellow puncturing tip, avoid damage to abdominal wall and suitable for second puncturing.



#### Available in 2 sizes

Code	Diameter	Length
T10100H	Ф10mm	100mm
T12100H	Ф12mm	100mm







### None-air inlet valve design

None-air inlet valve design, the overall structure is simple and the risk of air leakage is reduced.



#### Available in 3 sizes

Code	Diameter	Length
T05075N	Ф5mm	75mm
T05100N	Ф5mm	100mm
T05150N	Ф5mm	150mm

# Trocar Model-Single set



lmage	Model	Code	Diameter	Length
	Endopass™	T05075V T05100V T05150V	Ф5mm	75mm 100mm 150mm
a de la companya de l	Endopass™	T10100V T10150V	Ф10тт	100mm 150mm
o colors	Endopass™	T12100V T12150V	Ф12mm	100mm 150mm
and the second s	Quickpass™	T05075P T05100P	Ф5mm	75mm 100mm
and a second	Quickpass™	T10100P T12100P	Ф10mm Ф12mm	100mm
outres .	Easypass™	T05100S T10100S T12100S	Ф5mm Ф10mm Ф12mm	100mm
and the second s	Bluntpass™	T10100D T12100D	Ф10mm Ф12mm	100mm
0 000	Hassonpass™	T10100H T12100H	Ф10mm Ф12mm	100mm
	Without Valve ™	T05075V T05100V T05150V	Ф5mm	75mm 100mm 150mm

# **Packaging**





Model	Code	Obturator (Diameter)	Cannula (Unit)
	TV02	Ф5mm	Ф5mm x2
	TV06	Ф10mm	Ф10mm x2
Endopass™	TV08	Ф12mm	Ф12mm x2
Endopuss	TV28	Ф5mm Ф10mm Ф12mm	Ф5mm x1 Ф10mm x2 Ф12mm x2
	TV29	Ф5mm Ф10mm Ф12mm	Ф5mm x1 Ф10mm x3 Ф12mm x1
	TS02	Ф5mm	Ф5mm x1
	TS06	Ф10mm	Ф10mm x2
Easypass™	TS08	Ф12mm	Ф12mm <sub>х2</sub>
Lasypass	TS28	Ф5mm Ф10mm Ф12mm	Ф5mm x1 Ф10mm x2 Ф12mm x2
	TS29	Ф5mm Ф10mm Ф12mm	Ф5mm x1 Ф10mm x3 Ф12mm x1

# **Packaging**







# More Collection Available soon...

Balloon Trocar

Ф3mm Child Type

Ф15mm Trocar

**Thoracic Trocar** 

**Multipass Trocar** 

•

